

WIPE DISPENSER FOR FILLING WITH DISINFECTANTS

DESCO WIPES



DISPENSER BUCKET

1 / 4

Re-usable non-woven wipes dispenser **for filling
with DESCO WIPES non-woven wipes**

Flexible choice of disinfectant

**Simple
reprocessing
with Ultrazol
Oxy Sacht**




DESOTEX®
Wipes
System

WIPE DISPENSER FOR FILLING WITH DISINFECTANTS

DESCO WIPES



PRODUCT DESCRIPTION

2 / 4

DESCO WIPES is a reusable non-woven wipe dispenser system for disinfection and cleaning of surfaces of all types.

DESCO WIPES is made up of a reusable non-woven wipe dispenser system and corresponding rolls of non-woven wipes for refilling. It can be used in combination with both ready-to-use rapid disinfectants and ready-to-use application solution of disinfectant concentrates, making it very flexible.

Each refill pack comes with our ULTRASOL OXY SACHET to enable manual hygienic reprocessing before refilling the dispenser.

APPLICATIONS AND NOTES

Range of application

Refillable non-woven wipe dispenser system for impregnation, hygienic preparation and application of suitable surface disinfectants.

Application

For use with ready-to-use disinfectants or disinfectants prepared in accordance with manufacturer's instructions. We recommend the following products, which can be used for 28 days: OPTISAL PLUS, OPTISAL N, CLEANISEPT, BIGUANID FLÄCHE NR, DESCOSEPT SPEZIAL, DESCOSEPT PUR & ULTRASOL OXY. ¹

Wipe surfaces with the saturated DESCO WIPES until the surfaces are completely wetted and allow to dry for the entire contact time.

¹ The standing time of 28 days applies to all surface disinfectants of Dr. Schumacher GmbH as well as for the majority of all common products (follow manufacturer's instructions).

Users instruction

Filling: Remove the cover of the non-woven wipe dispenser and place with the outside on the work bench. Take a dry roll of non-woven wipes out of the packaging and insert into the non-woven wipe dispenser (alternatively: insert pre-saturated wipes with the bag). Pour the disinfectant solution slowly in a spiral over the fabric roll in the dispenser. Pull the first non-woven wipe from the middle of the roll of wipes and insert the front part of the wipe through the star-shaped dispensing system. Firmly close the non-woven wipe dispenser. When closing the dispensing opening, ensure that the wipe that has been drawn up is completely covered by the protective cover. Complete the enclosed cover label with the additional information and apply to the cover. After 30 minutes the wipes are ready to use. Discard the first wipe.

Each time the DESCO WIPES non-woven wipe dispenser is refilled, the system must undergo a preventive cleaning and disinfection. Along with possible chemical-thermal and thermal procedures, manual preparation using preparation wipes is possible. Each refill pack comes with our ULTRASOL OXY SACHET to enable manual hygienic reprocessing before refilling the dispenser. It is recommended to inspect the non-woven wipes dispenser before reprocessing. In case of wear marks cover/dispenser has to be replaced.

Cleaning: Thoroughly clean the insides and edge of the new or emptied dispenser with the cleaning wipe. Thoroughly clean the insides of the cover, the edges of the cover, inside surfaces and the edges of the dispenser opening. Thoroughly clean the outside of the cover and the outside surfaces of the dispenser.

Disinfection: Thoroughly wipe the insides and edge of the new or emptied dispenser with the disinfection wipe. Thoroughly wipe the insides of the cover, the edges of the cover, inside surfaces and the edges of the dispenser opening. Thoroughly wipe the outside of the cover and the outside surfaces of the dispenser.

Use disinfectants safely.

Always read label and product information before use.

WIPE DISPENSER FOR FILLING WITH DISINFECTANTS


DESCO WIPES




APPLICATIONS AND NOTES

3 / 4


Preparation

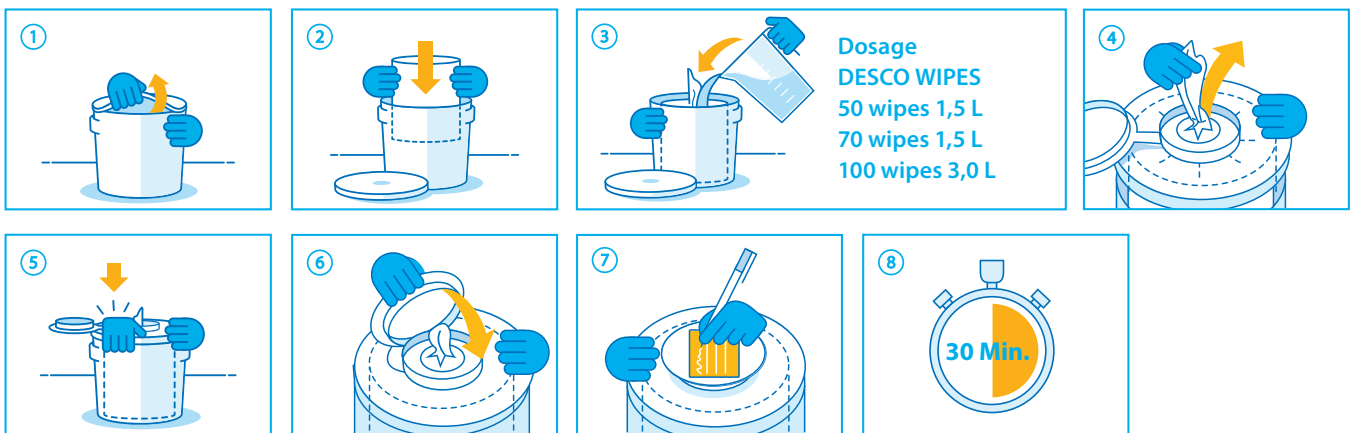

1 Cleaning

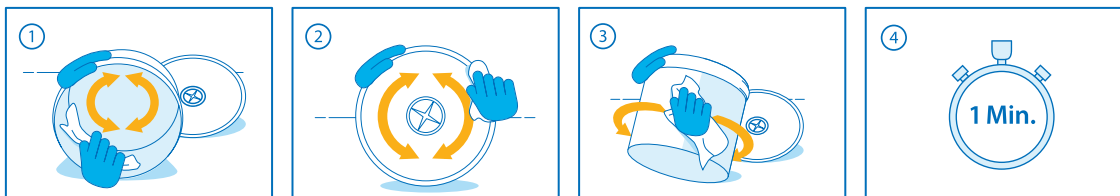

2 Disinfection

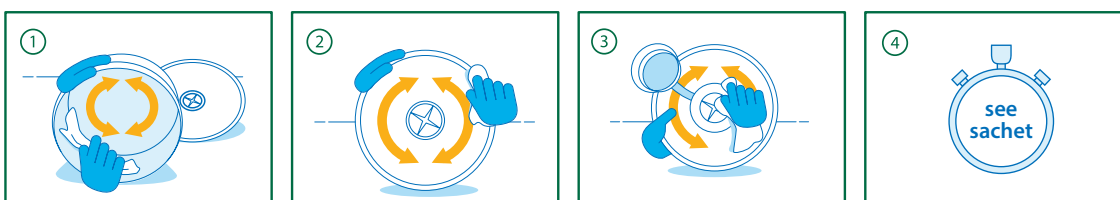

WIPE DISPENSER FOR FILLING WITH DISINFECTANTS

DESCO WIPES



PURCHASING INFORMATION

4 / 4

Product	Single unit	Delivery unit	Content	Size	REF
DESCO WIPES	dispenser bucket	1	empty	large	00-915-SD001
	dispenser bucket	6	empty	large	00-915-SD002
	dispenser bucket	1	empty	small	00-915-SKD001
	dispenser bucket	6	empty	small	00-915-SKD002
	wipes roll	3	100 wipes	30 x 32 cm	00-915-RD10003-01
	wipes roll	6	100 wipes	30 x 32 cm	00-915-RD10006-01
	wipes roll	3	70 wipes	25 x 25 cm	00-915-RD7003-01
	wipes roll	3	50 wipes	30 x 32 cm	00-915-RD5003-01
Application aids	wall holder	1		large	00-902-EIM
	wall holder	1		small	00-902-EIMK

CERTIFICATIONS



For an overview of our current certifications, scan the QR code to visit:

<https://www.schumacher-online.com/en/cert>

Dr. Schumacher is certified according to DIN EN 13485, DIN EN ISO 9001, DIN EN ISO 14001, BS OHSAS 18001, has a validated environment management system according to EMAS and is a member of IHO, VCI, BAH, DGSV and of the DGKH.