

NON-WOVEN WIPE DISPENSER SYSTEMS

ONE SYSTEM⁺
PLUS



SINGLE USE WIPE DISPENSER FOR FILLING WITH DISINFECTANT

1 / 4

Flexible choice of disinfectant

depending on the application area

Innovative design for easy handling

Completely **recyclable bucket**

Use after 15 minutes – **even in areas with a high risk of infection** – up to 60 days shelf life

- 18 % less plastic
- Folding carrying handle
- Optimized dispensing opening
- Compact design



NON-WOVEN WIPE DISPENSER SYSTEMS

ONE SYSTEM⁺ PLUS



PRODUCT DESCRIPTION

2 / 4

ONE SYSTEM PLUS is an innovative and proven disposable non-woven wipe dispenser system to be used as part of hygiene measures. Combined with suitable disinfectants, medical devices and medical equipment as well as surfaces of all types can be hygienically prepared.

ONE SYSTEM PLUS is quick and practical to use. The system is already fitted with a DESOTEX non-woven wipe roll so that users can carry out the required disinfection measure immediately after filling with the disinfectant and allowing to soak for 15 minutes. When used with Dr. Schumacher certain products,

ONE SYSTEM PLUS can be used over a period of up to 60 days instead of the 28 days usual for other products on the market. Once the non-woven wipe roll is finished, the system is not reprocessed but can be directly disposed of after any residue is removed. **ONE SYSTEM PLUS** is designed for single use and prevents possible contamination resulting from reprocessing and re-use. The system is therefore suitable for use in areas with an especially high risk of infection as well as the usual areas of application in the healthcare sector.

The innovative disposable non-woven wipe dispenser system is simple, safe and efficient. The construction and handling is similar to conventional non-woven wipe

dispensers but with a number of additional advantages: The dispenser opening is aligned with the non-woven wipe roll, giving the user quick access when removing the first wipe along with direct saturation of the entire roll and the resulting shortened preparation time. The filling process was designed to be intuitive and practical for the user. There are considerably fewer steps in the process for the user. The transparent design of the system means that the filling level can be checked at a glance, making it easy to plan the supply. The dispenser with its lightweight foldable construction is made of pure polypropylene and can be flattened for disposal without any special requirements.

The DESOTEX premium non-woven fabric in the system is a tear-resistant non-woven fabric and has been extensively reviewed by experts with the recommended disinfectants.

The **ONE SYSTEM PLUS** ensures optimal implementation of the desired disinfection measures with user-friendly handling, great cost effectiveness and demonstrated safe results.

Disinfectant	Shelf life after opening ¹	4-Field-Test ²	Use also in areas with special risk of infection ³
ULTRASOL OXY	60 days	•	•
OPTISAL PLUS	42 days	•	•
BIGUANID FLÄCHE NR	42 days	•	•
CLEANISEPT	42 days	•	•
OPTISAL N	42 days	•	•
DESCOCEPT PUR	42 days	•	•
DESCOCEPT SENSITIVE	42 days	•	•

¹ – In comparison with other conventional surface disinfectants with a shelf life after opening of 28 days and an impregnation time of 30 minutes, the tested surface disinfectants from Dr. Schumacher achieved a presaturation time of only 15 minutes and an increased shelf life after opening of the filled system of 42/60 days.

² – tested with all mandatory test organisms in line with EN 16615 (bactericidal/yeastcidal)

³ – in accordance with information from the Association for Applied Hygiene (VAH) Disinfectant Commission in co-operation with the "4+4 Working Group": On the use of wipe dispenser systems in areas with special infection risk, Hyg Med 2014; 39 – 9.

NON-WOVEN WIPE DISPENSER SYSTEMS

ONE SYSTEM⁺ PLUS



APPLICATIONS AND NOTES

3 / 4

Range of application

Disposable non-woven wipe dispenser system for impregnating, hygienically storing and using wipes with suitable surface disinfectants to disinfect and clean medical devices, medical equipment and surfaces of all types.

Application

For use with ready-to-use disinfectants or disinfectants prepared in accordance with manufacturer's instructions. We recommend the following products: OPTISAL PLUS, OPTISAL N, BIGUANID FLÄCHE NR, CLEANISEPT, DESCOSEPT SENSITIVE, DESCOSEPT PUR and ULTRASOL OXY. When handled correctly, a shelf life of up to 60 days can be achieved with the non-woven wipe dispenser system.

Unscrew the lid of the non-woven wipe dispenser.

Place the lid on a clean surface with the top facing down. Soak the roll of non-woven wipes already in the dispenser with 3.0 L of ready-to-use disinfectant solution. Open the protective cover of the lid. Pull the first non-woven wipe out from the middle of the roll of wipes. Draw the front part of the wipe through the star-shaped dispensing opening in the lid. Screw the lid back on and close the protective cover. When closing the dispensing opening, ensure that the wipe that has been drawn up is completely covered by the protective cover. Complete the enclosed label with the additional information. After 15 minutes saturation time, the non-woven wipes are ready to use. Discard the first wipe.

Wipe surfaces with the impregnated one system basic wipes until the surfaces are completely wetted and allow to dry for the entire contact time.

Single-use system! Dispose of any residual fluid, compact the dispenser with the lid open and place in the recycling system.

Application notes


Contents: The ONE SYSTEM PLUS contains 100 wipes 32 cm × 30 cm in size, made of tear-resistant non-woven material.

Disposal: The ONE SYSTEM PLUS disposable non-woven wipe dispenser is made of polypropylene and is designed for single use only. The ONE SYSTEM PLUS can be recycled completely after use. Discard any residual fluid before disposal.

According to the EU Medical Device Regulation, users/patients are obligated to report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the EU Member State in which the user/patient is established.

For professional use only by personnel with corresponding specialist knowledge according to national directives.

Read label before use.

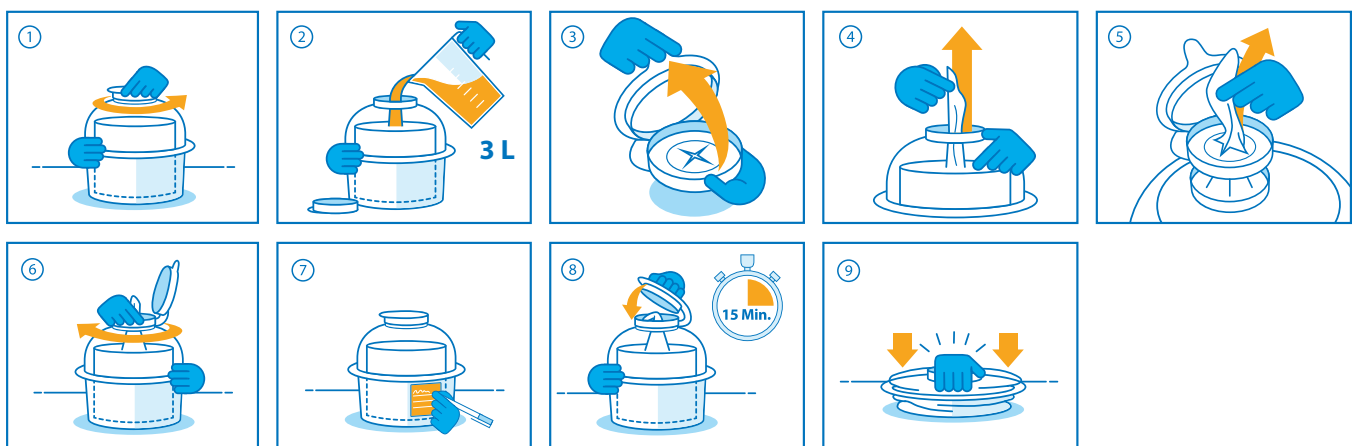

NON-WOVEN WIPE DISPENSER SYSTEMS

ONE SYSTEM⁺ PLUS



PURCHASING INFORMATION

4 / 4

Product	Single unit	Unit	Content	Wipe size	REF
ONE SYSTEM PLUS	Dispenser	2	100 Blatt	30 x 32 cm	00-915-OSEE100-01

The availability of the products and container sizes depend on a successful national registration.



RECOMMENDED DISINFECTANTS

When handled correctly, a shelf life of up to 60 days¹ can be achieved with the non-woven wipe dispenser system with the following products:

ULTRASOL OXY // OPTISAL PLUS // BIGUANID FLÄCHE NR // CLEANISEPT // OPTISAL N // DESCOSEPT PUR // DESCOSEPT SENSITIVE



Use disinfectants safely. Always read label and product information before use.

CERTIFICATIONS



For an overview of our current certifications, scan the QR code to visit:
<https://www.schumacher-online.com/en/cert>

Dr. Schumacher is certified according to DIN EN 13485, DIN EN ISO 9001, DIN EN ISO 14001, BS OHSAS 18001, has a validated environment management system according to EMAS and is a member of IHO, VCI, BAH, DGSV and of the DGKH.