

THERMO ALKA CLEAR



CLEAR-RINSING AGENT FOR THERMAL BEDPAN CLEANING MACHINES

1 / 3

Prevents calcification of bedpan washing machines


Promotes spotless drying of the bedpans

Low foaming



INSTRUMENTS / AUTOMATED REPROCESSING

THERMO ALKA CLEAR

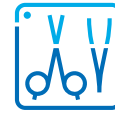

PRODUCT DESCRIPTION

2 / 3

THERMO ALKA CLEAR is suitable for alkaline thermal processing of bedpans in thermal washing machines for bedpans of all leading manufacturers. It ensures a spotless drying of the bedpan, reduces calcification inside the washing machine and is extremely lowfoaming.

RANGE OF APPLICATION

THERMO ALKA CLEAR is suitable for the alkaline thermal cleaning of bedpans in thermal cleaning machines for bedpans of all leading manufacturers. THERMO ALKA CLEAR reduces the calcification of the bedpan cleaning machines.

Application

Place materials into the machine according to instructions (e.g. hygiene plan). Every surface and cavity must be accessible for the washing liquor. Use THERMO ALKA CLEAR for post-rinsing according to the instructions of the bedpan manufacturer. Attention note for program setting: Do not add at temperatures below 30 °C.

According to the EU Medical Device Regulation, users/patients are obligated to report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the EU Member State in which the user/patient is established.

Application notes

Do not use with non-thermal and non-alkali resistant cleaning machines for bedpans.

Dosage

0.1 - 0.45 % (1 - 4,5 ml/L), depending on machine type and water hardness (see instructions of manufacturer).

Do not start addition below 30 °C. Dosage is to be done by qualified personnel.

Composition

15 - 30 % polyacrylate

Material compatibility

The use is not recommended for the following materials: Aluminum, copper, brass, Ethylene propylene diene (monomer) rubber (EPDM). We only recommend using under certain conditions for the following materials: Acrylonitril-butadiene-styrene (ABS).

THERMO ALKA CLEAR is suitable for the thermal treatment of Horcher, Arjo, ZK Hospital and Kodra bedpans.


For professional use only by personnel with corresponding specialist knowledge according to national directives.

The availability of the products and container sizes depend on a successful national registration.

Read label before use.

INSTRUMENTS / AUTOMATED REPROCESSING

THERMO ALKA CLEAR

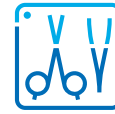

PURCHASING INFORMATION

3 / 3

Product	Single unit	Delivery unit	Content	REF
THERMO ALKA CLEAR	Flat canister	3	4,25 L	00-137-050



CERTIFICATIONS



For an overview of our current certifications, scan the QR code to visit:
<https://www.schumacher-online.com/en/cert>

Dr. Schumacher is certified according to DIN EN 13485, DIN EN ISO 9001, DIN EN ISO 14001, BS OHSAS 18001, has a validated environment management system according to EMAS and is a member of IHO, VCI, BAH, DGSV and of the DGKH.