

INSTRUMENTS / INSTRUMENT CARE

**SPEZIAL
ÖLSPRAY**




**OIL SPRAY BASED ON MEDICAL WHITE OILS FOR LUBRICATING MOVING PARTS
OF INSTRUMENTS AND MEDICAL DEVICES**

1 / 3



INSTRUMENTS / INSTRUMENT CARE

SPEZIAL ÖLSPRAY

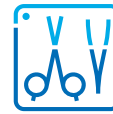

PRODUCT DESCRIPTION

2 / 3

SPEZIAL ÖLSPRAY based on pure medical white oils, is extremely compatible with common materials such as joints and taps on rigid endoscopes, as well as turbines, handpieces, contra-angles, and surgical instruments such as forceps and scissors. Furthermore the white oil complies with the DAB (German Pharmacopoeia) quality standard.

The special spray valve ensures a fine, even protective film. This determines the high vapor permeability of **SPEZIAL ÖLSPRAY** and allows safe sterilization in a steam sterilizer. **SPEZIAL ÖLSPRAY** increases the value retention of the materials through suitable lubrication of moving parts and effective corrosion protection.

APPLICATIONS AND NOTES

Application

For professional use only. Spray through turbines or handpieces and angle pieces twice daily or before each autoclaving process with **SPEZIAL ÖLSPRAY** for approx. 1 sec. and wipe off excess oil. Keep the can upright while spraying. Instruments: In a half-opened state, spray from approx. 25 cm distance. According to the EU Medical Device Regulation, users/patients are obligated to report any serious incident that has occurred in relation to the device to the manufacturer and the competent authority of the EU Member State in which the user/patient is established.

Application notes

Shake well before use. Always maintain this distance and spray on a thin oil film only.

Composition

Medical grade white oil (Paraffin perliquidum) according to DAB (German Pharmacopoeia), propellants propane/butane.

Material compatibility

If **SPEZIAL ÖLSPRAY** is used on latex or rubber these materials may swell.

Precautionary and hazard statements

Extremely flammable aerosol. Pressurised container: May burst if heated. Keep away from heat, hot surfaces, sparks, open flames and other ignition sources. No smoking. Do not spray on an open flame or other ignition source. Do not pierce or burn, even after use. Protect from sunlight. Do not expose to temperatures exceeding 50 °C/ 122 °F. Dispose of contents/ container to an approved waste disposal plant.

Read label before use.

INSTRUMENTS / INSTRUMENT CARE

SPEZIAL
ÖLSPRAY



PURCHASING INFORMATION

3 / 3

Product	Single unit	Delivery unit	Content	REF
SPEZIAL ÖLSPRAY	300 ml aerosol can	12	300 ml	00-108-003

National information may differ. For further information, please contact our subsidiary or your local dealer.
The availability of the products and container sizes depend on a completed national registration.



CERTIFICATION



For an overview of our current certifications, scan the QR code to visit:
<https://www.schumacher-online.com/en/cert>

Dr. Schumacher is certified according to DIN EN 13485, DIN EN ISO 9001, DIN EN ISO 14001, BS OHSAS 18001,
has a validated environment management system according to EMAS and is a member of IHO, VCI, BAH, DGSV and of the DGKH.